

Wound and Skin Care Solution with Stabilizing Preservatives

Life-Altering Advance in Wound Care

Fast Acting

- Cleans infections in one step
- Safe to use around eyes, nose, mouth and open wounds
- Accelerates healing of minor wounds and skin irritation
- Non-toxic and non-irritating



Wound and Skin Care for Those Who are Active



Puracyn Wound and Skin Care is optimal for immediate care of cuts, abrasions and lacerations. Always have Puracyn on hand so you are prepared when someone in your family gets injured. Puracyn cleans infections, accelerates healing in one step while it is also non-toxic and non-irritating.

Add Puracyn to your family's hygiene and wound care program. It is excellent for cleaning injuries including cuts, lacerations, skin abrasions and skin irritations. Outside or inside the house, Puracyn, based on revolutionary Microcyn[®] Technology, should be your first line of protection. It does not contain antibiotics, steroids or alcohol and no prescription is required. Safe to use around eyes, ears, nose and mouth!

99.999% Inactivation in-solution in just 30 Seconds


- | | |
|--------------------------------|--------------------------|
| ✗ MRSA - Staphylococcus aureus | ✗ Escherichia coli |
| ✗ VRE - Enterococcus faecalis | ✗ Salmonella |
| ✗ Acinetobacter baumannii | ✗ Pseudomonas aeruginosa |
| ✗ Candida albicans | |



To Order:

Call 1-866.318.3116 or visit www.puracyn.com

The Science of Microcyn® Technology



Does not contain Antibiotics, Steroids or Alcohol

Puracyn® OTC


Kills 99.99% of Staph, MRSA, Strep and Other Infectious Germs in 30 Seconds In-Solution

NEW!

Wound and Skin Care
Germ-Killing Solution with Stabilizing Preservatives

FAST ACTING!

- Cleans infections in one step
- Safe to use around eyes, nose, mouth and open wounds
- Accelerates healing of minor wounds and burns
- Non-toxic and non-irritating



INDICATIONS: Puracyn Wound and Skin Care is optimal for immediate care (early treatment) of cuts, abrasions, lacerations and minor burns.

DIRECTIONS FOR USE: Apply the Puracyn Wound and Skin Care directly on the wound to flush away dirt and debris. Continue application until the wound is thoroughly cleaned. If dressing is used, saturate it with Puracyn and apply directly over the wound. Dressing size should extend beyond the wound margins. Flush the wound with Puracyn and apply a fresh, moistened dressing 1-5 times daily as needed. Soak adherent dressings with Puracyn to facilitate removal.

WARNING: Not for injection. In case of deep or puncture wounds, animal bites or serious burns, consult a health care professional. If condition worsens or persists for longer than seven days, consult a health care professional. May cause a temporary stinging sensation. Keep out of reach of children.

HANDLING: Non-irritating, non-toxic and non-sensitizing to the skin and eyes. No special handling precautions required. Biodegradable and environmentally safe solution.

STORAGE: Store at room temperature away from direct sunlight and heat. **DO NOT FREEZE.** Do not use after the expiration date.

STABILITY: Product has a shelf life up to 18 months from date of manufacture when stored at room temperature.


DISPOSAL: No special disposal considerations are required.

INGREDIENTS: Electrolyzed Water (99.97%), Sodium Chloride (NaCl) 0.023%, and preservatives: Sodium Hypochlorite (NaOCl) 0.004%, and Hypochlorous Acid (HOCl) 0.003%.

Puracyn® OTC
Manufactured by Oculis Innovative Sciences for Innovacyn, Inc.
3546 N. Riverside Ave, Rialto, CA 92377
www.puracyn.com • 1-866-318-3116

Made in the USA

Lot Number: _____
Expiration date: _____



Part Number: 6002; 80001 - Pump v1.0

8 fl. oz. (236 mL)

Microcyn® Technology Mode of Action

Microcyn Technology is the only safe and stable wound cleanser containing hypochlorous acid, a compound naturally produced by the human immune system (neutrophils). The production of hypochlorous acid is generated via a proprietary manufacturing process that enables the solution to remain stable. Clinical studies to date have demonstrated the safety of the product in cleaning and debriding a wide range of wounds. See FDA-cleared label indications for Puracyn.

Performance Testing

Puracyn Wound and Skin Care Solution with stabilizing preservatives has been evaluated by the USP<51> (in-solution) antimicrobial effectiveness test and in time-kill studies. It has demonstrated effectiveness by reducing the microorganisms in-solution (listed to right) in-solution and passes the USP antimicrobial effectiveness test for a category one product.

Name of Organism	Log Reduction (30 sec.)	Time to Kill	Percent Reduction
Staphylococcus aureus MRSA	6.34	30 seconds	99.9999%
Enterococcus faecalis VRE	6.36	30 seconds	99.9999%
Staphylococcus aureus	6.23	30 seconds	99.9999%
Escherichia coli	5.70	30 seconds	99.9997%
Acinetobacter baumannii	6.37	30 seconds	99.9999%
Bacteroides fragilis	7.64	30 seconds	99.9999%
Candida albicans	6.33	30 seconds	99.9999%
Enterobacter aerogenes	6.09	30 seconds	99.9999%
Enterococcus faecium VRE - MDR	6.51	30 seconds	99.9999%
Haemophilus influenzae	5.18	30 seconds	99.9993%
Klebsiella oxytoca MDR	6.05	30 seconds	99.9999%
Klebsiella pneumoniae	6.14	30 seconds	99.9999%
Micrococcus luteus	5.84	30 seconds	99.9999%
Proteus mirabilis	6.20	30 seconds	99.9999%
Pseudomonas aeruginosa	5.81	30 seconds	99.9998%
Serratia marcescens	6.00	30 seconds	99.9999%
Staphylococcus epidermidis	6.02	30 seconds	99.9998%
Staphylococcus haemolyticus	5.91	30 seconds	99.9999%
Staphylococcus hominis	5.45	30 seconds	99.9996%
Staphylococcus saprophyticus	5.96	30 seconds	99.9999%
Streptococcus pyogenes	6.72	30 seconds	99.9999%

Puracyn Wound and Skin Care with germ-killing preservatives is not FDA cleared as an in-wound antimicrobial. Ingredients, warnings and instructions for use are identified in adjacent label.



To Order:
Call 1-866.318.3116 or visit www.puracyn.com


Innovacyn, Inc.
3546 N. Riverside Ave., Rialto, CA 92377
Part No. 85005 v1.0 MKT31088-02 Rev. 2

Puracyn[®] OTC

Wound and Skin Care Solution with Stabilizing Preservatives

Life-Altering Advance in Wound Care

YOUR FIRST LINE OF DEFENSE

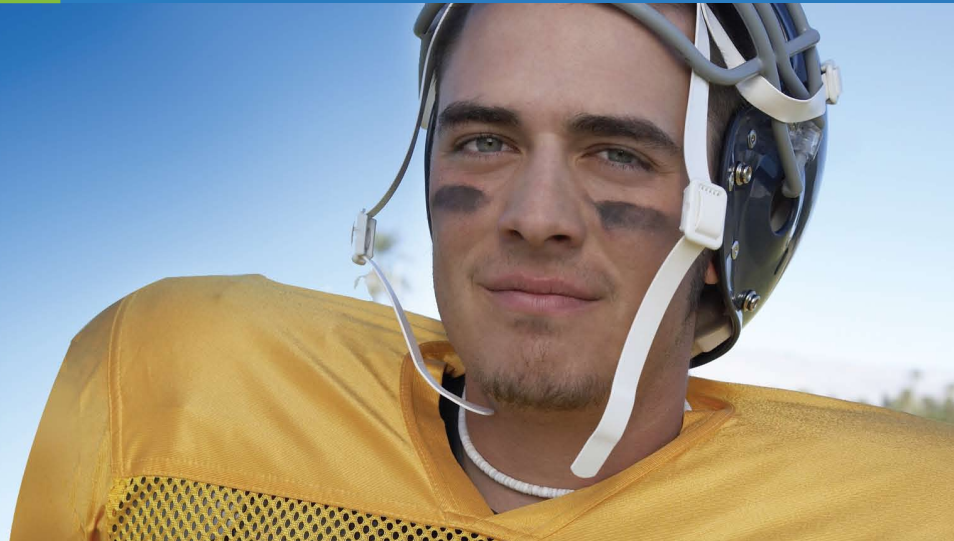

Wound and Skin Care for Athletes



Athletes who participate in contact sports are particularly susceptible to skin and soft-tissue damage caused by turf burns, abrasions, shaving and chafing. These injuries can be further complicated by MRSA infections, which can be extremely dangerous and in some cases, deadly. All contact sports are the most commonly affected sports, due to the greater amount of skin-on-skin contact and the large number of athletes congregating in common locker rooms.

Add Puracyn to your team's hygiene and wound care program. Excellent for cleaning sports injuries including lacerations, skin abrasions and minor skin irritations. On the field and off, Puracyn, should be your first line of defense. Puracyn does not contain antibiotics, steroids or alcohol and no prescription is required. Safe to use around the eyes, nose and mouth!

99.999% Inactivation in-solution in just 30 seconds

- ✕ MRSA - Staphylococcus aureus
- ✕ Escherichia coli
- ✕ VRE - Enterococcus faecalis
- ✕ Salmonella
- ✕ Acinetobacter baumannii
- ✕ Pseudomonas aeruginosa
- ✕ Candida albicans



To Order:

Call 1-866.318.3116 or visit www.puracyn.com

NEW!

The Science of Microcyn® Technology

Does not contain Antibiotics, Steroids or Alcohol

Puracyn OTC

Kills 99.99% of Staph, MRSA, Strep and Other Infectious Germs in 30 Seconds In-Solution

Puracyn OTC
Manufactured by Oculus Innovative Sciences for Innovacyn, Inc.
3546 N. Riverside Ave, Rialto, CA 92377
www.puracyn.com • 1-866-318-3116

Made in the USA

Lot Number:
Expiration date:



Part Number: 6002; 80001 - Pump v1.0

Puracyn OTC

Wound and Skin Care
Germ-Killing Solution with Stabilizing Preservatives

FAST ACTING!

- Cleans infections in one step
- Safe to use around eyes, nose, mouth and open wounds
- Accelerates healing of minor wounds and burns
- Non-toxic and non-irritating



8 fl. oz. (236 mL)

INDICATIONS: Puracyn Wound and Skin Care is optimal for immediate care (early treatment) of cuts, abrasions, lacerations and minor burns.

DIRECTIONS FOR USE: Apply the Puracyn Wound and Skin Care directly on the wound to flush away dirt and debris. Continue application until the wound is thoroughly cleaned. If dressing is used, saturate it with Puracyn and apply directly over the wound. Dressing size should extend beyond the wound margins. Flush the wound with Puracyn and apply a fresh, moistened dressing 1-5 times daily as needed. Soak adherent dressings with Puracyn to facilitate removal.

WARNING: Not for injection. In case of deep or puncture wounds, animal bites or serious burns, consult a health care professional. If condition worsens or persists for longer than seven days, consult a health care professional. May cause a temporary stinging sensation. Keep out of reach of children.

HANDLING: Non-irritating, non-toxic and non-sensitizing to the skin and eyes. No special handling precautions required. Biodegradable and environmentally safe solution.

STORAGE: Store at room temperature away from direct sunlight and heat. **DO NOT FREEZE.** Do not use after the expiration date.

STABILITY: Product has a shelf life up to 18 months from date of manufacture when stored at room temperature.

DISPOSAL: No special disposal considerations are required.

INGREDIENTS: Electrolyzed Water (99.97%), Sodium Chloride (NaCl) 0.023%, and preservatives: Sodium Hypochlorite (NaOCl) 0.004%, and Hypochlorous Acid (HOCl) 0.003%.

Microcyn® Technology Mode of Action

Microcyn Technology is the only safe and stable wound cleanser containing hypochlorous acid, a compound naturally produced by the human immune system (neutrophils). The production of hypochlorous acid is generated via a proprietary manufacturing process that enables the solution to remain stable. Clinical studies to date have demonstrated the safety of the product in cleaning and debriding a wide range of wounds. See FDA-cleared label indications for Puracyn.

Performance Testing

Puracyn Wound and Skin Care Solution with stabilizing preservatives has been evaluated by the USP<51> (in-solution) antimicrobial effectiveness test and in time-kill studies. It has demonstrated effectiveness by reducing the microorganisms in-solution (listed to right) in-solution and passes the USP antimicrobial effectiveness test for a category one product.

Name of Organism	Log Reduction (30 sec.)	Time to Kill	Percent Reduction
Staphylococcus aureus MRSA	6.34	30 seconds	99.9999%
Enterococcus faecalis VRE	6.36	30 seconds	99.9999%
Staphylococcus aureus	6.23	30 seconds	99.9999%
Escherichia coli	5.70	30 seconds	99.9997%
Acinetobacter baumannii	6.37	30 seconds	99.9999%
Bacteroides fragilis	7.64	30 seconds	99.9999%
Candida albicans	6.33	30 seconds	99.9999%
Enterobacter aerogenes	6.09	30 seconds	99.9999%
Enterococcus faecium VRE - MDR	6.51	30 seconds	99.9999%
Haemophilus influenzae	5.18	30 seconds	99.9993%
Klebsiella oxytoca MDR	6.05	30 seconds	99.9999%
Klebsiella pneumoniae	6.14	30 seconds	99.9999%
Micrococcus luteus	5.84	30 seconds	99.9999%
Proteus mirabilis	6.20	30 seconds	99.9999%
Pseudomonas aeruginosa	5.81	30 seconds	99.9998%
Serratia marcescens	6.00	30 seconds	99.9999%
Staphylococcus epidermidis	6.02	30 seconds	99.9998%
Staphylococcus haemolyticus	5.91	30 seconds	99.9999%
Staphylococcus hominis	5.45	30 seconds	99.9996%
Staphylococcus saprophyticus	5.96	30 seconds	99.9999%
Streptococcus pyogenes	6.72	30 seconds	99.9999%

Puracyn Wound and Skin Care with germ-killing preservatives is not FDA cleared as an in-wound antimicrobial. Ingredients, warnings and instructions for use are identified in adjacent label.



To Order:
Call 1-866.318.3116 or visit www.puracyn.com

Innovacyn, Inc.
3546 N. Riverside Ave., Rialto, CA 92377
Part No. 85000 v1.2 MKT31088-03 Rev. 2